



# Table Top Sterilizer (Front Loading)



## Table Top Sterilizer (Front Loading)



**OBROMAX**

**ACMAS Technologies Inc.**

CE An ISO 9001:2000 Company

Sales Office: 312-313, Vardhman Capital Mall, L.S.C. 10, Gulabi Bagh, Delhi-110052, INDIA  
Tel: +91-11-23654603 (M) +91-9313971681, 9311039044, 9555351619 Fax: +91-11-23654603  
Email: [acmastechnologies@gmail.com](mailto:acmastechnologies@gmail.com) | [obromax.delhi@gmail.com](mailto:obromax.delhi@gmail.com)  
[www.acmastechnologies.com](http://www.acmastechnologies.com) | [www.acmas.in](http://www.acmas.in)

**ACMAS Technologies Inc.**

CE An ISO 9001:2000 Company

**OBROMAX**


# Table Top Sterilizer (Front Loading)



- This is fully automatic front loading autoclave with dry cycle.
- Electronic Circuit inside the autoclave with programmable integrated circuit chip.
- With dry cycle after sterilization without opening the door.
- With Digital Display of all the operations in sequence and timer and temperature.
- The front control panel has a touch type key switches and also have LEDs to indicate the operation going inside the autoclave.
- The autoclave is programmed to perform all necessary operations required in sterilization.
- With additional solenoids valve provided to shunt out the air from the chamber at the initial stage.
- The internal chamber have two number of trays made of stainless steel having size 12" x 8"



# Table Top Sterilizer (Front Loading)

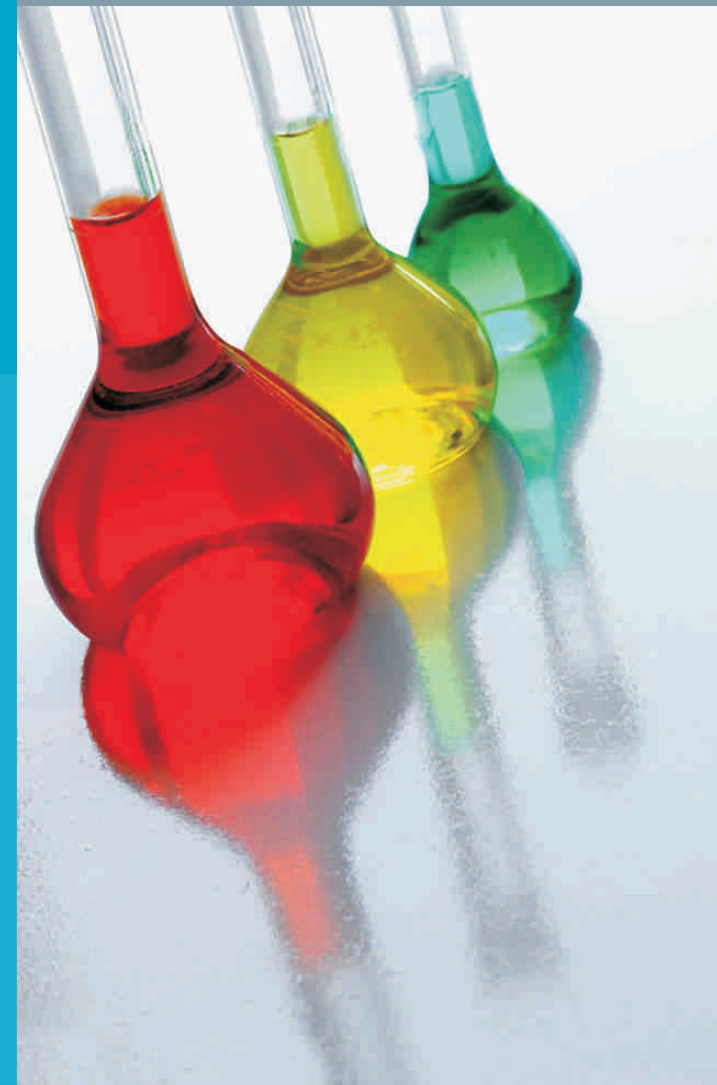

- Sealing Ring made of silicone rubber.
- One sterilization cycle.
- Sterilization with dry cycle completes with single touching of start button.

**Overall Size :**

**Length :**

480mm, Width : 470mm, Height : 380mm

**Capacity :**

12 Litre

**Power Consumption :**

1 KW

**Gross Weight :**

24 KG